



TDE Series ultra-low temperature freezers



TDE Series ultra-low temperature freezers

Thermo Scientific™ TDE Series –80°C ultra-low temperature (ULT) freezers, powered by the H-Drive system, are designed for efficiency. As a result, the series delivers a foundational ULT experience with the reliability and performance needed to protect your work—and designed to help meet your sustainability goals.

Designed to keep your samples safe

- Worry-free storage; with low peak variation, samples are protected throughout the entire cabinet¹
- Fast door opening recovery; as little as 2 minutes for chest freezers
- Secure freezer settings with sample security
- General-purpose models for laboratory use, and FDA Class II Medical Device models for preserving blood and blood products by storing them at cold or freezing temperatures
- Backed by our Global Service and Support team

Limit environmental impact and lower HVAC costs

- Energy use as low as 4.4 kWh/day, depending on size and temperature
- Less heat emission lowers HVAC costs
- SNAP and EU F-gas regulation compliant, natural hydrocarbon refrigerants
- Environmentally friendly, polyurethane foam with non-HFC blowing agent is compliant with SNAP and EU F-gas regulations
- Manufactured in an award-winning, zero-waste-to-landfill facility

User-friendly and quiet

- H-Drive Information Center (HIC) touchscreen user interface
- Quiet enough to reside directly in the laboratory; at just 50 dBA, the sound is comparable to a home refrigerator³
- 5-year parts and labor standard warranty*

* Restrictions and exclusions apply. Warranty is void if the unit is not maintained as set forth in the operation and service manuals.

1. Typical freezer data based on internal testing with freezer setpoint at –80°C and ambient temperature at 20°C. Performance data vary across freezer size and electrical specification.

2. FDA Class II Models are listed under product code KSE, regulation number 864.9700.

3. All data based on internal data unless otherwise noted. Data on file.



A rich heritage of innovation

From our first ULT freezer in 1940 to today, the combined experience of REVCO, Forma Scientific, Jewett, Harris, and other families of products continues to drive our ability to deliver industry-leading innovations that help accelerate science and health care.

Easy-to-access controls with setpoint security


- Easy to use: plug in and go, no setup required
- HIC combines easy-to-read temperature with simple capacitive touch buttons
- Advanced security features: alarm settings, setpoint security passcode, and padlock-compatible doors



Heart—the heart is the health status for the freezer

Easy-to-read alerts

- Alarm bell**—indicates visual and audible alarms
- Door**—illuminates during a door ajar alarm
- Wrench**—generic service warning
- Thermometer**—indicates when the cabinet temperature exceeds either warm alarm or cold alarm setpoints

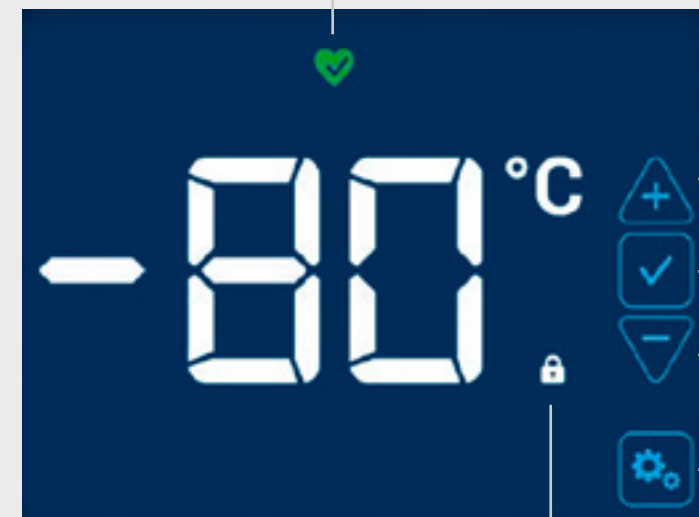

Easy-to-navigate controls

- Plus**—increases the value of the selected setting
- Check mark**—saves a change to the selected value
- Minus**—decreases the value of the selected setting
- Settings**—represents the various settings

Setpoint security—protect your settings with a 3-digit numeric code

Remote monitoring

All TDE Series ULT freezers can be managed with optional Thermo Scientific™ Smart-Vue Pro™ remote monitoring. Smart-Vue Pro enables you to view the status of your equipment—and samples—24/7 through a cloud or premise-based portal. Smart-Vue Pro also provides 21 CFR Part 11 validated data logging, eliminating a need for additional resources to perform manual processes.

Smart-Vue Pro Benefits

- Monitoring-Sensors collect temperature, relative humidity, differential pressure, door openings and more, giving you greater situational awareness of your equipment
- Traceability-Auditable, secure system logs with complete, history of actions and events to eliminate manual logging
- Sample security-programmable parameters with alerts and notifications can let you know if freezer contents are at risk
- Installation, IQ/OQ, calibration and maintenance plans available

For more information on Smart-Vue Pro, please visit thermofisher.com/smartvuepro



Verified sustainability

TDE Series ULT freezers have received the ACT Environmental Impact Label from the nonprofit organization, My Green Lab. The label addresses the need of both scientists and procurement specialists for clear, third party-verified information about the environmental impact of laboratory products. The criteria for the ACT label were developed with input from industry experts and external stakeholders, and the label process entails third-party verification of the sustainable impacts of a product, its operation, and its end of life. For additional information on the ACT label, go to thermofisher.com/actlabel.

User-friendly cabinet design

- Ergonomic design: single-hand operation with lockable door handle for easy and secure access to samples
- Flat top workspace for quick access to hands-free operation
- Casters for easy mobility; all are lockable to safely secure unit into place
- Easy-to-access USB for simplified service



A sound environment

You know that the constant noise created by compressors can compromise communications and create a less-than-ideal working environment. TDE Series ULT freezers are up to 10x quieter than our previous generation of ULT freezers.*



Decibel scale of common sounds

Noise source	Decibel level
Freeway traffic, vacuum cleaner	70
Conversation in restaurant, office, background music	60
TDE Series ULT freezers	50
Library	40
Quiet rural area	30
Rustling leaves	22
Breathing	10

* Thermo Scientific™ TDE60086FA model compared to Thermo Scientific™ Revco™ UXF60086FA model.

TDE Series -80°C ULT upright freezers

Thermo Scientific™ TDE Series -80°C ULT upright freezers, powered by the H-Drive system, are designed for efficiency. As a result, the series enables a foundational ULT experience with the reliability and performance needed to protect your work—and designed to help meet your sustainability goals.



Performance data

Model series	Energy consumption* (kWh/day)	Peak variation from -80°C setpoint (°C)**	60 second door opening recovery time (minutes)	Warm-up time (-80°C to -50°C) (minutes)†	Noise level (dBA)
TDE30086FV	9.2	+6.3 / -3.6	11	268	51.5
TDE40086FV	10.4	+5.7 / -3.5	15	271	50.8
TDE50086FV	10.3	+4.2 / -3.8	16	312	50
TDE60086FV	11.9	+5.8 / -3.3	21	352	50.5

* Energy consumption: Typical data based on internal testing with freezer setpoint at -80°C and 20°C ambient temperature. ** Peak variation from setpoint: Typical data based on internal testing with freezer setpoint at -80°C and 20°C ambient temperature. † Warm-up time: Typical freezer data based on internal testing with empty freezer setpoint at -80°C and ambient temperature at 20°C.

Upright freezers (-50°C to -86°C)

General-purpose model	Capacity		Electrical	Amps/breaker (plug type)	Max shelf weight lbs. (kg)	Interior dimensions H x W x D in. (cm)	Exterior dimensions H x W* x D in. (cm)	Shipping weight lbs. (kg)
	Cu. ft. (L)	Cryoboxes (2 in.)						
TDE30086FV	14.9 (422)	300	230V/50Hz	5.9/10 (European)	125 (57)	51.2 x 17.8 x 28.3 (130.1 x 45.2 x 71.9)	78.0 x 23.0 x 38.5 (198.1 x 58.4 x 97.7)	712 (323)
TDE40086FV	19.4 (549)	400	230V/50Hz	5.3/10 (European)	165 (75)	51.2 x 23.1 x 28.3 (130.1 x 58.8 x 71.9)	78.0 x 28.3 x 38.5 (198.1 x 71.9 x 97.7)	687 (312)
TDE50086FV	24.1 (682)	500	230V/50Hz	5.9/10 (European)	205 (93)	51.2 x 28.8 x 28.3 (130.1 x 73.0 x 71.9)	78.0 x 34.0 x 38.5** (198.1 x 86.4 x 97.7)	734 (333)
TDE60086FV	28.8 (816)	600	230V/50Hz	6.1/10 (European)	245 (111)	51.2 x 34.4 x 28.3 (130.1 x 87.3 x 71.9)	78.0 x 39.6 x 38.5** (198.1 x 100.6 x 97.7)	820 (372)

* Allowance for handle and hinge, add 4.2 in. (10.7 cm) to cabinet width. ** Doorway width clearance is 34.5 in. (87.6 cm) for TDE50086 and TDE60086 models. For 5 inner door version, refer to "Options & Accessories" table on page 9

TDE Series -80°C ULT chest freezers

Thermo Scientific™ TDE Series -80°C ULT chest freezers, powered by the H-Drive system, enable the same excellent experience as our ULT upright freezers, with an additional benefit of enhanced performance. The combination of our ULT chest freezers' rapid door open recovery and tight peak variation make them perfect for frequent entry or long-term storage applications.



Performance data

Model Series	Energy consumption* (kWh/day)	Peak Variation from -80°C setpoint ©**	60-second door opening recovery time (minutes)	Warm-up time (-80°C to -50°C) (minutes)†	Noise level (dBA)
TDEC06386FV	4.9	+4.5/-2.1	2	156	50.2
TDEC25286FV	8.5	6.7/-3.1	2	154	51.6
TDEC39686FV	10.8	5.9/-3.5	2	156	53.7

Specifications for chest freezers (-50°C to -86°C)

General-purpose model	Capacity		Electrical	Amps/breaker (plug type)	Interior dimensions H x W x D in. (cm)	Exterior dimensions H x W x D in. (cm)	Shipping weight lbs. (kg)
	Cu. ft. (L)	Cryoboxes (2 in.)					
TDEC06386FV	3 (85)	63	230V/50Hz	5.6/16 (European)	16.4 X 18.4 x 18.4 (41.7 x 46.7 x 46.7)	43.7 x 28.8 x 32.7 (111.0 x 73.2 x 83.1)	380 (176)
TDEC25286FV	13 (360)	252	230V/50Hz	5.2/16 (European)	27.8 x 42.4 x 18.5 (70.6 x 107.7 x 47.0)	40.5 x 72.3 x 33.1 (102.9 x 183.6 x 84.1)	710 (322)
TDEC39686FV	20 (566)	396	230V/50Hz	5.5/16 (European)	27.8 x 66.4 x 18.5 (70.6 x 168.7 x 47.0)	40.5 x 96.3 x 33.1 (102.9 x 244.6 x 84.1)	860 (390)

TDE Series -40°C ULT upright freezers

Thermo Scientific™ TDE Series -40°C ULT freezers, powered by the H-Drive system, bring the excellent experience you'd expect from our -80°C freezers but tuned specifically for the performance required between -10°C and -40°C.



Performance data

Model Series	Energy consumption* (kWh/day)	Peak Variation from -80°C setpoint ©**	60-second door opening recovery time (minutes)	Warm-up time (-80°C to -50°C) (minutes)†	Noise level (dBA)
TDE40040FV	4.4	+3.2 / -1.6	19	309	50.8
TDE50040FV	5.3	+1.8 / -2.2	22	328	50.7
TDE60040FV	6.2	+1.7 / -1.9	29	339	48.1

* Energy consumption: Typical data based on internal testing with freezer setpoints at -40°C and 20°C ambient temperature.
 ** Peak variation from setpoint: Typical data based on internal testing with freezer setpoint at -40°C and 20°C ambient temperature.
 † Door opening recovery: Typical data—internal recovery testing is -38.5°C. During the door opening, all inner and outer doors are opened for 60 seconds.
 ‡ Warm-up time: Typical data based on internal testing with empty freezer setpoint at -40°C and ambient temperature at 20°C.
 For specific model performance data, please refer to the individual technical specification sheet for the desired model.

Specifications for upright freezers (-10°C to -40°C)

General-purpose model	Capacity		Electrical	Amps/breaker (plug type)	Max shelf weight lbs. (kg)	Interior dimensions H x W x D in. (cm)	Exterior dimensions H x W* x D in. (cm)	Shipping weight lbs. (kg)
	Cu. ft. (L)	Cryoboxes (2 in.)						
TDE40040FV	19.4 (549)	400	230V/50Hz	3.2/10 (European)	165 (75)	51.2 x 23.1 x 28.3 (130.1 x 58.8 x 71.9)	78.0 x 28.3 x 38.5 (198.1 x 71.9 x 97.7)	642 (292)
TDE50040FV	24.1 (682)	500	230V/50Hz	3.4/10 (European)	205 (93)	51.2 x 28.8 x 28.3 (130.1 x 73.0 x 71.9)	78.0 x 34.0 x 38.5** (198.1 x 86.4 x 97.7)	668 (304)
TDE60040FV	28.8 (816)	600	230V/50Hz	3.4/10 (European)	254 (111)	51.2 x 34.4 x 28.3 (130.1 x 87.3 x 71.9)	78.0 x 39.6 x 38.5** (198.1 x 100.6 x 97.7)	765 (348)

TDE Series ULT upright freezers— options and accessories

Options (field installations must be performed by a qualified professional)

Description	Installation	Model series			
		TDE30086	TDE40086	TDE50086	TDE60086
LN₂ Backup System Maintains temperature down to -80°C with liquid nitrogen. Models for TDE30086 are stand-alone.	Factory	N/A	LN4567	LN4567	LN4567
	Field	FLN300D* FLN300V** FLN300W†	FLN456LED* FLN456LEV**	FLN456LED* FLN456LEV**	FLN456LED* FLN456LEV**
CO₂ Backup System Maintains temperature down to -67°C with CO ₂ . Models for TDE30086 are stand-alone.	Factory	N/A	CO4567	CO4567	CO4567
	Field	FCO300D* FCO300V** FCO300W†	FCO456LED* FCO456LEV**	FCO456LED* FCO456LEV**	FCO456LED* FCO456LEV**
TSX Pro Series User Interface Option (only available for upright models)	Factory	UI34567H	UI34567H	UI34567H	UI34567H
Chart Recorder (Inkless) Seven-day inkless recorder Range: -115°C to +50°C (+5°C resolution) Dia.: 6 in. (15.2 cm)	Factory	N/A	CR4567	CR4567	CR4567
	Field	N/A	FCR400LE	FCR500LE	FCR600LE
Chart Recorder (Ink) Seven-day ink recorder Range: -100°C to +38°C (+2°C resolution) Dia.: 6 in. (15.2 cm) Models for TDE30086 and TDE40086 are stand-alone.	Factory	N/A	CR4567	CRP4567	CRP4567
	Field	6383-6D* 6383-6V** 6383-6W†	6383-6D* 6383-6V** 6383-6W†	FCRP500LE	FCRP600LE
Revco Five Inner Door Configuration	Factory	5D300	5D400	5D500	5D600
Stainless Steel Freezer Interior	Factory	SS34567	SS34567	SS34567	SS34567
Twist Lock Plug NEMA L6-15P for 208-230V/60Hz	Factory	US230V15ATL	US230V15ATL	US230V15ATL	US230V15ATL
	Field	TL230V15A	TL230V15A	TL230V15A	TL230V15A
Twist Lock Plug NEMA L6-20P for 208-230V/60Hz	Factory	US230V20ATL	US230V20ATL	US230V20ATL	US230V20ATL
	Field	TL230V20A	TL230V20A	TL230V20A	TL230V20A

* D = 208-230V/60Hz (NEMA 6-15 plug). ** V = 230V/50Hz (European plug). † W = 230V/50Hz (UK plug).

Accessories

Description	Model series			
	TDE30086	TDE40086	TDE50086	TDE60086
Air Filter (Replacement)	AF34567LE	AF34567LE	AF34567LE	AF34567LE
Alarm Delay Module Eliminates nuisance alarms due to intermittent or transitory conditions. (The alarm conditions must exist for a user-defined period before the signal is broadcast.)	6903	6903	6903	6903
Battery Backup (Replacement)	400159	400159	400159	400159
Chart Paper, for Inkless Recorder, 50/Pack	N/A	N/A	6185	6185
Chart Paper, for Ink Recorder, 50/Pack	17020	17020	17020	17020
Cryo Gloves, Medium	4425	4425	4425	4425
Cryo Gloves, Large	4426	4426	4426	4426
Seismic Restraint Kits	TF-ULT300	TF-ULT400	TF-ULT500	TF-ULT600
Shelf Kits: Include One Shelf and Clips	SK300	SK400	SK500	SK600







Cord description	Model
NEMA 5-20 125V/20A	US125V20A
NEMA 6-20 230V/15A	US230V15A
Argentina 230V/16A	AR230V16A
Australia 230V/16A	AU230V16A
Brazil 230V/16A	BR230V16A
China 230V/16A	CH230V16A
Denmark 230V/16A	DK230V16A
Continental EU 230V/16A	EU230V16A
India 230V/16A	IN230V16A
Israel 230V/16A	IS230V16A
Italy 230V/16A	IT230V16A
Japan 125V/20A	JP125V20A
Japan 250V/15A	JP250V15A
Switzerland 230V/16A	SW230V16A
UK 230V/13A	UK230V13A

TDE Series ULT chest freezers— options and accessories

Accessories

Description	Installation	Model series		
		TDEC06386	TDEC25286	TDEC39686
LN ₂ Backup System, maintains temperature down to -80°C with liquid nitrogen	Factory	332353G01	332353G05	332353G05
CO ₂ Backup System, maintains temperature down to -70°C with CO ₂	Factory	332353G02	332353G06	332353G06
Chart Recorder (Inkless), 7-day inkless recorder Range: -115°C to +50°C (+5°C resolution) Diameter: 6 in. (15.2 cm)	Factory	332791G01	332166G01	332166G01
Chart Recorder (Ink), 7-day ink recorder Range: -115°C to +50°C (+5°C resolution) Diameter: 6 in. (15.2 cm)	Factory	332791G02	332166G02	332166G02



Description	Model No.
Replacement Air Filter Chest Freezers (quantity of 5)	332297H02
Replacement Backup Battery	400159TF
Chart Paper, for Inkless Recorder, 50/pack	90520H07
Chart Paper, for Ink Recorder, 50/pack	17020
Cryo Gloves, medium	4425
Cryo Gloves, large	4426

	Description/dimensions H x W x D in. (cm)	Freezer capacity cu. ft. (liters)	Boxes/rack (plates)	Racks/freezer (plates)	Boxes/freezer (plates)	Cat. No.
	2 Box Rack with locking rod, 5.5 x 15.2 x 5.5 (14 x 38.6 x 14)	3 (85)	7	9	63	398184
	3 Box Rack, 5.5 x 5.5 x 16 (14 x 14 x 40.6)	3 (85)	5	9	45	398185
	2 Box Rack with locking rod, 5.5 x 5.6 x 26.5 (14 x 14.2 x 67.3)	13 (360)	12	21	252	398186
		20 (566)	12	33	396	
	3 Box Rack, 9.4 x 5.5 x 26.5 (14 x 14 x 67.3)	13 (360)	8	21	168	398187
		20 (566)	8	33	264	
	27 Position Microplate Rack, 5.3 x 3.5 x 25.8 (13.5 x 8.9 x 65.5)	13 (360)	(27)	(33)	(891)	398188
		20 (566)	(27)	(51)	(1377)	
	Wire Basket with Rollers* 11.8 x 16.5 x 12.3 (29.9 x 41.9 x 31.2)	13 (360)	NA	5	NA	114048
		20 (566)	NA	9	NA	



*Basket rail set 189000 is required.

Racks for boxes and microplates

Racks for 5 cm boxes

Description	Dimensions H x W x D cm	Storage capacities	Model series				Cat. No.
			TDE30086	TDE40086	TDE50086	TDE60086	
Sliding drawer rack	30.2 x 14.0 x 68.3	Boxes/rack	25	25	25	25	 920090
		Racks/shelf	3	4	5	6	
		Racks/freezer	12	16	20	24	
		Boxes/freezer	300	400	500	600	
Adjustable side access rack	29.5 x 13.7 x 67.9	Boxes/rack	25	25	25	25	 1950520
		Racks/shelf	3	4	5	6	
		Racks/freezer	12	16	20	24	
		Boxes/freezer	300	400	500	600	

Racks for 7.62 cm boxes




Description	Dimensions H x W x D cm	Storage capacities	Model series				Cat. No.
			TDE30086	TDE40086	TDE50086	TDE60086	
Sliding drawer rack	30.2 x 14.0 x 68.3	Boxes/rack	15	15	15	15	 920091
		Racks/shelf	3	4	5	6	
		Racks/freezer	12	16	20	24	
		Boxes/freezer	180	240	300	360	
Adjustable side access rack for 3 in. boxes	29.5 x 13.7 x 67.9	Boxes/rack	15	15	15	15	 1950521
		Racks/shelf	3	4	5	6	
		Racks/freezer	12	16	20	24	
		Boxes/freezer	180	240	300	360	

Racking shelf kits

Description	Storage capacities	Model series			
		TDE30086	TDE40086	TDE50086	TDE60086
	Cat. No.	RSK300SD4	RSK400SD4	RSK500SD4	RSK600SD4
Sliding drawer racks with boxes for one shelf (12.7 x 12.7 x 5.0 cm boxes with 100-count cell dividers)	Racks included	3	4	5	6
	Boxes included	75	100	125	150

Services

Racks for standard or deep-well microplates

Description	Dimensions H x W x D cm	Storage capacities	Model series				Cat. No.
			TDE30086	TDE40086	TDE50086	TDE60086	
Side access rack with locking rod	30.2 x 14.0 x 65.3	Plates/rack	147	147	147	147	 1950592
		Racks/shelf	3	4	5	6	
		Racks/freezer	12	16	20	24	
		Plates/freezer	1,764	2,352	2,940	3,528	
Side access microplate rack	30.2 x 14.0 x 65.3	Plates/rack	105	105	105	105	 1950523
		Racks/shelf	3	4	5	6	
		Racks/freezer	12	16	20	24	
		Plates/freezer	1,260	1,680	2,100	2,520	
Sliding drawer microplate rack	30.2 x 14.0 x 66.3	Plates/rack	35	35	35	35	 1950642
		Racks/shelf	3	4	5	6	
		Racks/freezer	12	16	20	24	
		Plates/freezer	420	560	700	840	

Boxes and dividers

Description	No. of dividers	Quantity	Color	Cat. No.
Thermo Scientific™ Nalgene™ polycarbonate 5.08 cm box for 1.2 to 2.0 mL vials	81	24/case	White	5026-0909
5.08 cm fiberboard box	0	12/pack	White	5954
5.08 cm fiberboard box	81	1	White	820002
5.08 cm fiberboard box	100	1	White	820109
7.62 cm fiberboard box	0	12/pack	White	5956
7.62 cm fiberboard box	81	1	White	820003
5.08 cm fiberboard box	0	12/pack	Multiple	1950878
5.08 cm fiberboard box	0	12/pack	Pink	1950879
5.08 cm fiberboard box	0	12/pack	Red	1950880
5.08 cm fiberboard box	0	12/pack	Yellow	1950881
5.08 cm fiberboard box	0	12/pack	Green	1950882
5.08 cm fiberboard box	0	12/pack	Blue	1950883
5.08 cm fiberboard box	0	12/pack	Black	1950884

Fiberboard dividers

Description	Quantity/pack	Cat. No.
Grid dividers, hold 100 x 12 mm tubes/vials	12	5958
Grid dividers, hold 81 x 13 mm tubes/vials	12	6212
Grid dividers, hold 64 x 14 mm tubes/vials	12	5960
Grid dividers, hold 49 x 16 mm tubes/vials	12	5959

Certificate of calibration

Control, alarm, and temperature recording systems are systematically calibrated at the factory for accurate performance. An optional certificate of calibration to an NIST (National Institute of Standards and Technology, formerly National Bureau of Standards) traceable instrument includes a model number and serial number certificate (single-point calibration). Multiple-point calibration is also available. Calibration certificates must be requested when ordering. For details, contact your sales representative.

Validation testing program

Validation testing services ordered in advance of shipment may accelerate in-house installation qualification (IQ) and operation qualification (OQ) programs. Validation services include documentation and performance reports, compressor and subsystem component specifications, microprocessor flow diagrams, software specifications, and more. For more details, contact your sales representative.

Services description

Temperature mapping of 24 hours at setpoint to include:

- Mapping data graphed at each location
- Thermocouple placement diagram
- Test equipment list
- Average uniformity and stability at specified cycle
- Minimum, average, and maximum temperatures at each thermocouple location
- Uniformity and stability at each location

74554TA

Temperature mapping of 24 hours (traceable instruments and raw data) at setpoint to include:

- Mapping data graphed at each location
- Thermocouple placement diagram
- Test equipment list utilizing Kaye Validator™ thermal validation equipment
- Calibration certificates of test equipment
- Thermocouple calibration reports
- Average uniformity and stability at specified cycle
- Minimum, average, and maximum temperatures at each thermocouple location
- Uniformity and stability at each location

7456TA

Certificate of verification* 6916-2

Factory quality control test log report* 6916-4

Factory quality control test temperature graphs* 6916-5

Calibration certificate** 260029

Certificate of conformance** 260043

Unit-specific calibration certificate** 260049

*For upright models

**For chest models, to be ordered separately



Also available: TSX Series ULT freezers

Designed to support sustainability and meet the highest protection standards for storage at ultra-low temperatures.

Find out more at thermofisher.com/cold

thermo scientific